



No Lifting Required!


Patient Transfer Board

Roller-Free Roller Board - Better For Patients, Easier On Nurses, and Inexpensive

- **No Rollers** - Less traumatic transfer for patients
- **Comfortable** - Soft but sturdy board is easy on patients
- **Lightens the Load** - No lifting required, even with your heaviest patients
- **Lightweight** - Very easy to use
- **No Special Training Requirements** - Uses same or similar transfer technique
- **Radiolucent** - You can image with the board under the patient
- **Fowler-Position Transfer** - Long board allows for transfer in seated position
- **Foldable** - Long board folds for easy storage and set-up
- **Straps** - Carry this lightweight device on your shoulder and keep your hands free



Makes Fowler-Position Transfer a Snap!

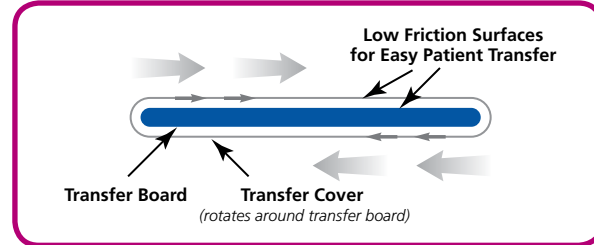

Ultra-Lightweight!
Short Board is 3 lbs. (1.5 kg)
Long Board is 6.3 lbs. (3 kg)

#A-83000 Long Patient Transfer Board
71" x 20" (180 cm x 50 cm)

#A-83001 Short Patient Transfer Board
39" x 16" (100 cm x 40 cm)


#A-83002 Disposable Covers (Case of 50)
(sized for Long Boards)



#A-83001



#A-83000

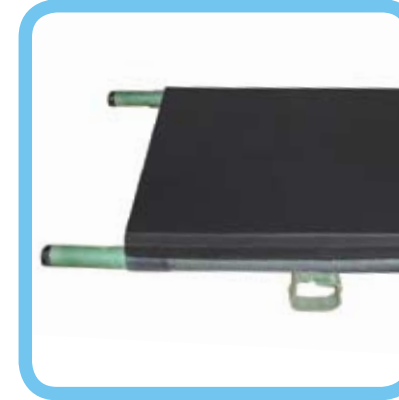

Warrior Evacuation Litter Pad

The Warrior Evacuation Litter Pad or "WELP" was developed to meet the needs of injured Service Members returning to the U.S on extended Aero-Medical Evacuation flights. Engineered specifically to meet the high standards set by the U.S Air Force's "Air Mobility Command" the WELP was approved by the Air Force for aviation use but is easily adaptable for a variety of applications.

The WELP features a durable non-skid bottom material that prevents sliding of the pad over the litter surface during steep ascent and rapid deceleration frequently experienced in tactical take off and landing. Baffled vents on the edge of the pad are provided to diminish pad swelling in the event of in-flight emergency rapid decompression.

Innovative features and manufacturing techniques such as ultrasonically welded seams and specialized sealants reduce the risk of the interior foams wicking blood and body fluids during use. Even during prolonged use, WELP's layered foam system coupled with a specialized patient interface surface aid in reducing the risk of soft tissue injuries such as Decubitis and Neuropathy.

Specs:

Color: Black

Dimensions: 72"x 23"x 2"

Weight: 10.5 lbs

* Not intended as a flotation device.

Ordering information:

Allen Medical Systems Catalog Number- H-223-48-223-A1

National Stock Number - 6530-01-548-0262

* Call for contract pricing.



#H-223-48-223-A1 Warrior Evacuation Litter Pad